

Azlon®

REUSABLE PLASTIC LABWARE

Bridging the
language gap

Look inside
for more!

Azlon® Multilingual Wash Bottles



Azlon® Multilingual Wash Bottle Range

The Azlon® Multilingual Wash Bottle range displays chemical information in four European languages and also meets the latest requirements in European safety legislation – plus all of the regular features and benefits of this popular laboratory essential.



Chemical names in English, French, German and Spanish eliminate the need for multiple wash bottles across European language speaking countries.

Chemical formulae and CAS number for easy identification of solvent contained.

Globally Harmonised System reflects the latest changes in European legislation including new pictograms.

Signal word (i.e. 'Danger'), to identify the severity of the hazard statements, replace the old risk phrases.


Key Features

- Wide mouth for easy filling
- Colour-coded lids and contents
- Precision jet tip for fine control
- Choice of indelibly printed solvent labels
- Custom printing available on request



Precision jet tip for fine control.

Wide mouth for easy filling and colour-coded cap and bottle for easy recognition of solvent.

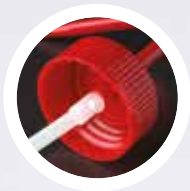


Available with **Driplok**® vapour venting technology featuring a unique, patented two-way valve.

Safety information in four languages – English, French, German and Spanish.

Solving the Solvent Problem

Volatile solvents can create a build-up of gas, forcing solvent up and out of the spout. This causes a steady drip of highly dangerous, flammable liquid.

The **Driplok**® venting wash bottle has a unique, patented two-way valve. It allows the vapour to escape when pressure builds up and stops air getting into the liquid.

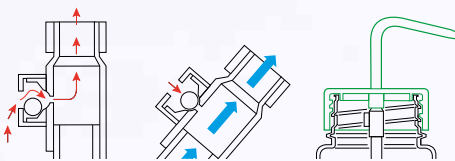

The Result:

No more hazardous drips or dribbles and an accurate, continuous dispensing stream.



Azlon® Round Integral Wash Bottle

We are proud to introduce the new Round Integral Wash Bottle to the Azlon® range. With this design, you can be sure to dispense liquids in a safe and efficient manner with extra protection from contamination.

The new Azlon® Round Integral Wash Bottle displays chemical names and key safety information in four European languages and also meets the latest requirements in European safety legislation – plus all of the regular features and benefits of this popular laboratory essential.



Traditional wash bottles consist of several components – a bottle, a cap and an inner tube. The inner tube can be a source of contamination when refilling, as it is typically placed on a lab bench surface.

Chemical names in English, French, German and Spanish eliminate the need for multiple wash bottles across European language speaking countries.

Chemical formulae and CAS number for easy identification of solvent contained.

Globally Harmonised Symbols reflects the latest changes in European legislation including new pictograms.

Signal word (i.e. 'Danger'), to identify the severity of the hazard statements, replaces the old risk phrases.



A range of designs to meet your needs



Vented cap with PTFE membrane allows venting and prevents pressure build-ups.

Precision jet tip ensures fine control.

Moulded graduations for approximate volume measurement.

Integral inner tube prevents bench-derived contamination when refilling.

Azlon® DIY and Tracker Wash Bottle

The Azlon® DIY and Tracker Wash Bottles are supplied complete with a unique write-on, wipe-off labelling facility.

Easy to reuse time and time again, the Tracker Wash Bottle is specially designed for full traceability. The bottles even meet the stringent labelling requirements of the Clinical Pathology Accreditation Body.



- Allows full traceability of solutions
- Tracker Wash Bottle meets the labelling requirements of the CPA
- All information is clearly visible to all users, at all times
- Easy to reuse
- Driplok® vapour venting valve ensures safe use with volatile solvents
- Self-adhesive COSHH labels supplied with all bottles



The information is easily applied using a solvent based pen and can be removed by simply wiping the bottle clean with acetone. **No more do-it-yourself labels!**

Technical and Ordering Information

Bottles, Wash, Multilingual, Safety Venting, Wide Neck, LDPE

- Multilingual Wide Neck Wash Bottles feature the patented **DripLok®** vapour venting valve that prevents solvent drips
- Indelibly printed in up to four languages (please refer to the **Region Key*** for a list of the languages that are printed on each wash bottle)
- Choose from five solvent labels all including the new Globally Harmonised Symbols, NFPA diamond, risk phrases and safety information
- Fitted with colour-coded, one-piece PP closures that maintain their shape in use
- Closures feature precision jet-tip aperture for fine control



Product Code	Capacity (ml)	Print Label	Closure Colour	Country Region	Pack Quantity
WGW531VTML	250	Acetone	Red	Region 1	5
WGW535VTML	250	Ethanol	Orange	Region 1	5
WGW534VTML	250	I.M.S	Yellow	Region 1	5
WGW533VTML	250	Isopropanol	Blue	Region 1	5
WGW532VTML	250	Methanol	Green	Region 1	5
WGW537VTML	500	Acetone	Red	Region 1	5
WGW541VTML	500	Ethanol	Orange	Region 1	5
WGW540VTML	500	I.M.S	Yellow	Region 1	5
WGW539VTML	500	Isopropanol	Blue	Region 1	5
WGW538VTML	500	Methanol	Green	Region 1	5
WGW550VTML	500	I.D.A	Yellow	Region 1	5
WGW537VTMLS	500	Acetone	Red	Region 2	5
WGW541VTMLS	500	Ethanol	Orange	Region 2	5
WGW539VTMLS	500	Isopropanol	Blue	Region 2	5
WGW538VTMLS	500	Methanol	Green	Region 2	5
WGW537VTMLR	500	Acetone	Red	Region 3	5
WGW541VTMLR	500	Ethanol	Orange	Region 3	5
WGW539VTMLR	500	Isopropanol	Blue	Region 3	5
WGW538VTMLR	500	Methanol	Green	Region 3	5
WGW537VTMLRE	500	Acetone	Red	Region 4	5
WGW541VTMLRE	500	Ethanol	Orange	Region 4	5
WGW539VTMLRE	500	Isopropanol	Blue	Region 4	5
WGW538VTMLRE	500	Methanol	Green	Region 4	5

*Region Key For Multilingual Wash Bottles

Region 1				Region 3			
English	French	German	Spanish	English	Russian	Polish	German
Region 2				Region 4			
English	Norwegian	Danish	Swedish	English	Arabic	French	

Azlon® Multilingual Round Integral Wash Bottle

- Multilingual Wide Neck Wash Bottle
- 500ml capacity
- Indelibly printed in English, French, German and Spanish
- Choose from five solvent labels all including the new Globally Harmonised Symbols, NFPA diamond, risk phrases and safety information
- Fitted with colour-coded, one-piece PP closure that maintains its shape in use
- Designed with integral spout to prevent contamination



Product Code	Print Label	Closure Colour	Pack Quantity
WGW801P	Distilled Water	White	5
WGW802P	Acetone	Red	5
WGW803P	Methanol	Green	5
WGW804P	Isopropanol	Blue	5
WGW805P	Ethanol	Orange	5

Bottles, Wash, DIY, Wide Neck, LDPE

- Do-it-yourself **DripLok®** Wash Bottles with blank write-on panel
- Indelibly printed in three languages (English, French and German)
- Supplied with adhesive hazard labels



Product Code	Description	Card Type	Pack Quantity
WGW001VTML	250ml Custom/DIY with white cap	Vapour-venting	5
WGW002VTML	500ml Custom/DIY with white cap	Vapour-venting	5
WGW020VTML	1000ml Custom/DIY with white cap	Vapour-venting	5
WGWPEN01	Indelible black marker pen for Custom/DIY wash bottles	–	2
WGW004VTML	500ml Custom/DIY with mixed colour caps	Vapour-venting	5
WGW001PML	250ml Custom/DIY with white cap	Non-venting	5
WGW002PML	500ml Custom/DIY with white cap	Non-venting	5
WGW020PML	1000ml Custom/DIY with white cap	Non-venting	5

Bottles, Wash, Tracker, LDPE

- Wash bottle designed with complete labelling facility to ensure full traceability
- Meets industry standards including those required by the Clinical Pathology Accreditation body
- Ideal for use in the clinical and general hospital market
- Unique write-on, wipe-off panel
- **DripLok®** vapour venting cap version available – for ultimate safety with volatile solvents
- Supplied with a set of COSHH labels

i Recordable information includes: bottle contents, first use and expiry dates, lot number, COSHH information, originator and date of receipt/preparation.



Product Code	Capacity (ml)	Cap Type	Closure Colour	Pack Quantity
WTR250P	250	Non-venting	Green	5
WTR250VT	250	Vapour-venting	Red	5
WTR500P	500	Non-venting	Green	5
WTR500VT	500	Vapour-venting	Red	5



SciLabware is a company of **DWK Life Sciences**

Email us with your questions at enquiries@scilabware.com

SciLabware Limited, Unit 4, Riverside 2,
Campbell Road, Stoke-on-Trent, ST4 4RJ, UK

Tel: +44 (0)1782 444406 | Fax: +44 (0)1782 940436